

**Product Code**

UTS-0170 Mechanical End-Over-End Shaker  
 UTS-0171 Gas Jar with a Rubber Cover for UTS-0170  
 UTS-0172 Ground Glass for UTS-0171

The UTS-0170 Mechanical End-Over-End Shaker is used for the determination of the particle density by the gas jar method and the particle size distribution by sedimentation.

End-Over-End Shaker is used to rotate two gas jars with rubber covers at about 50 rpm.


The 1 L capacity gas jar is made of plexiglass.

UTS-0171 Gas Jar with a rubber cover and a ground glass should be ordered separately.

Models for 220-240V 50 Hz, 1 ph.	UTS-0170-T
Models for 110-120V 60 Hz, 1 ph.	UTS-0170-
Models for 220-240V 60 Hz, 1 ph.	N UTS-0170-K

**Standards**

BS 1377:2

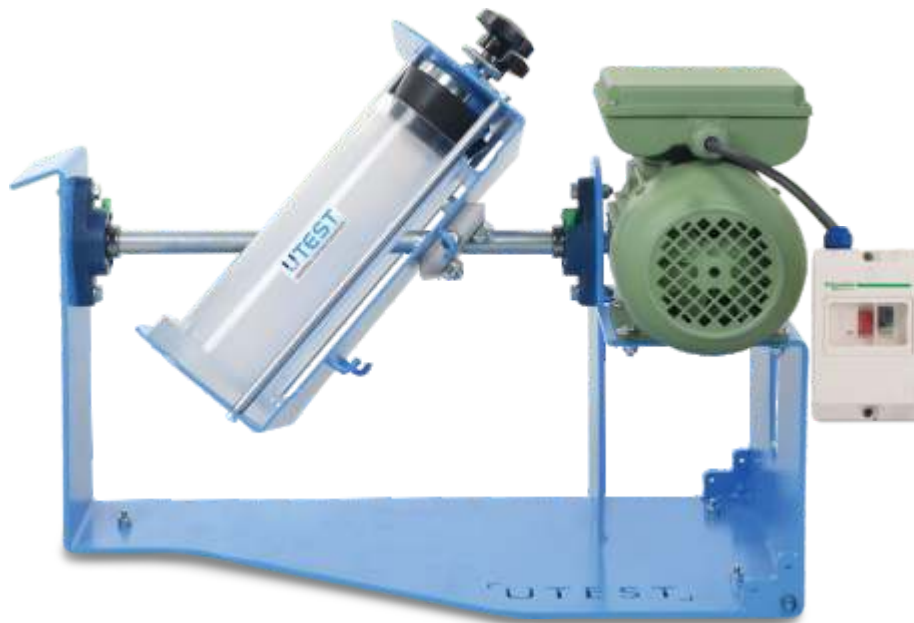

UTS-0170



UTS-0172

UTS-0171

Product Code	Dimensions	Weight (approx.)	Power
UTS-0170	500x1030x450	33 kg	180 W
UTS-0171	mm	0,60	
	120x120x280	kg	
	mm		